solutions™

Dialite® Adjusting, Finishing & Polishing System



225.00.638

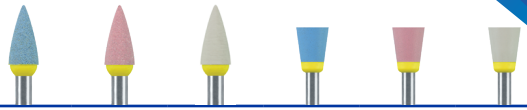


Figure No.	W16DGHP	W16DMHP	W16DHP	W17DGHP	W17DMHP	W17DHP
Head Size ($^1/_{10}$ mm)	040	040	040	060	060	060
Head Length (mm)	10.0	10.0	10.0	7.5	7.5	7.5



Figure No.	8801LDF	8369DF	8856DF
Head Size ($^1/_{10}$ mm)	014	025	016
Head Length (mm)	–	5.5	8.0


Order online at BrasselerUSA.com
US: 800.841.4522 or fax 888.610.1937
Canada: 800.363.3838 or fax 866.330.4454

BRASSELER USA®

©2021 Brasseler USA. All rights reserved.

B-3744B-REV3-5022180L0-02.21

solutions™ 225.00.638 Dialite® Adjusting, Finishing & Polishing System

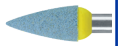	Order No.	Head Size	Grit	Usage
	W16DCHP.RA	040	C	Polisher - Occlusal Surface
	W16DMHP.RA	040	M	Polisher - Occlusal Surface
	W16DHP.RA	040	F	Polisher - Occlusal Surface
	W17DCHP.RA	060	C	Polisher - Labial/Lingual/Buccal

	Order No.	Head Size	Grit	Usage
	W17DMHP.RA	060	M	Polisher - Labial/Lingual/Buccal
	W17DHP.RA	060	F	Polisher - Labial/Lingual/Buccal
	8801LDF.FG	014	F	Diamond - Endo Access/ Small Adjustments
	8369DF.FG	025	F	Diamond - Occlusal Surface
	8856DF.FG	016	F	Diamond - Contact Surfaces/ Crown Removal

TIPS & TECHNIQUES: Polishers produce best results when using all three steps in order (coarse, medium and fine); Use Copious amounts of water when using the diamonds.

Complimentary items: LD Grinders for heatless grinding (LD11C/M)

Power Suggestion: Use a Brasseler Forza F1 Electric Attachment at low speed for polishers and Brasseler Forza F5 high speed attachment for diamonds.

Sterilization: Autoclave only

Recommended Speeds: 5,000-7,000 rpms

BRASSELER USA®