



BRASELER
USA®



EVERclean
MAINTENANCE MACHINE



BY YOUR SIDE

EVERclean

Designed with convenience in mind, the EverClean Maintenance Machine is the perfect solution for today's busy practice. Combined lubrication and purging cycles lead to a quick and easy, hassle-free maintenance that keeps your handpieces running like new. Affordable, fast, and thorough, the EverClean is the ideal maintenance solution for all your handpiece needs.


EverClean™ Maintenance Machine
5028208U0

· 2-Year Warranty




Complete with three **Ports** for air or electric attachments. Adapters available for most handpieces and attachments.

Adapters sold separately

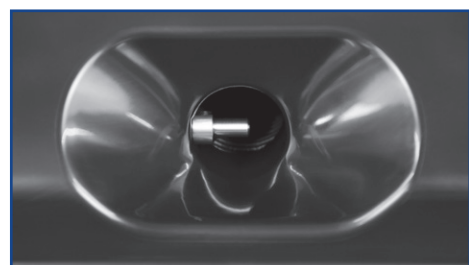

Internal Suction Valve creates vacuum effect to capture larger oil droplets from the cleaning chamber, making clean up even easier.

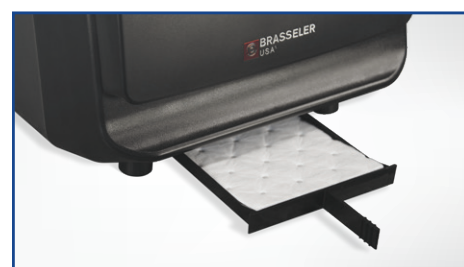

Removable **Oil Tray** with replaceable **Oil Sheets** collect excess lubrication and debris during the purging cycle.



Simple **Display Panel** allows you to easily choose desired handpiece stations, activate the chuck cleaning cycle, and adjust cycle time as needed.



Dedicated **Chuck Cleaning Station** to improve the performance and further extend the life of your handpiece.



EverClean™ Spray
5028210U0

· 500ml spray can
· Over 400 uses per can
Spray can sold separately

Removable **Door Shield** catches excess oil and debris for easy cleaning and can be replaced when needed.



EverClean™ Spray (500 ml)
5028210U0



EverClean™ Oil Sheets (10/Pk)
5028209U0



EverClean™ Door Shield (3/Pk)
5028220U0

HIGH SPEED AIR ADAPTERS



EverClean™ Kavo® Adapter
5028211U0



EverClean™ W&H® Adapter
5028212U0



EverClean™ Brasseler/NSK® Adapter
5028213U0



EverClean™ Sirona® Adapter
5028216U0



EverClean™ Midwest® Adapter
5028217U0



EverClean™ Star® Adapter
5028218U0

LOW SPEED & ELECTRIC ADAPTERS



EverClean™ E-Type Adapter
5028214U0



EverClean™ Midwest® Adapter
5028219U0

KaVo® is a registered trademark of Kaltenbach & Voigt GmbH, Germany. Midwest® is a registered trademark of DENTSPLY International. NSK® is a registered trademark of Nakanishi, Inc. Sirona® is a registered trademark of Dentsply Sirona, Inc. Star® is a registered trademark of DentalEZ, Inc. W&H® is a registered trademark of W&H Dentalwerk Bürmoos GmbH, Austria.

Visit our website at BrasselerUSA.com

To order call 800.841.4522 or fax 888.610.1937. In Canada call 800.363.3838

