




K0408

Evanesce™

Universal Direct Composite Contouring, Finishing and Polishing System

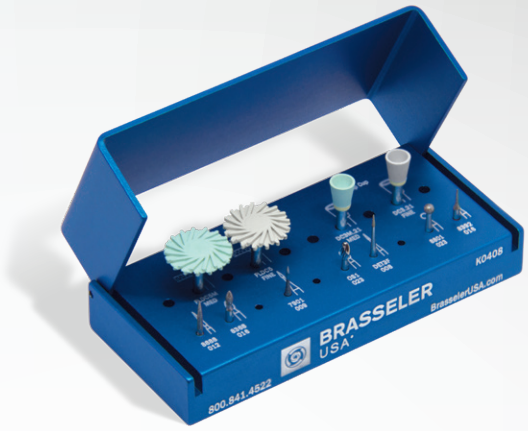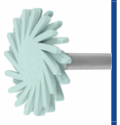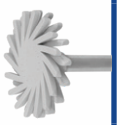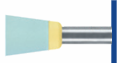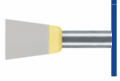

	Product Name	US No.	Order No.	Head Size ($^1/_{10}$ mm)	Cutting Length (mm)	Grit/ Blade	Usage
	Peter Brasseler Series™ Flame Diamond	–	8888	012	8.0	F	Create anterior surface texture Form developmental depressions Bevel Class IV tooth preparations
	Peter Brasseler Series™ Football Diamond	–	8368	016	3.5	F	Create surface texture; perikymata/imbrication lines Remove lingual flash Finish lingual/occlusal surface Occlusal adjustment
	Flame Trimming & Finishing Carbide	7901	H246	009	3.6	12 Blade	Create surface texture Form developmental depressions, finish lines Finish gingival margins Create occlusal surface anatomy Remove flash and overhang
	OS1 Football Trimming & Finishing Carbide	OS1 (7408)	H379	023	4.2	12 Blade	Occlusal adjustment Remove excess restorative material Finish lingual/occlusal surface Create secondary anatomy
	Peter Brasseler Series ™ Needle Diamond	DET3F	132F	008	3.0	F	Refine gingival margin
	Peter Brasseler Series™ Round Diamond	–	8801	023	–	F	Remove bubbles Occlusal adjustment Finish lingual/occlusal surface
	Peter Brasseler Series ™ Needle/Mosquito Diamond	–	8392	016	5.0	F	Finish gingival/incisal embrasures Remove flash and overhang
	Diacomp™ Feather Lite™ Multi-Use Premium 2-Step Polisher	–	FLDC5M	140	1.6	M	Pre-polish all surfaces
	Diacomp™ Feather Lite™ Multi-Use Premium 2-Step Polisher	–	FLDC5	140	1.6	F	Final polish all surfaces
	Diacomp™ Cup Multi-Use Premium 2-Step Polisher	–	DC3M.21	060	7.5	M	Pre-polish all surfaces
	Diacomp™ Cup Multi-Use Premium 2-Step Polisher	–	DC3.21	060	7.5	F	Final polish all surfaces

Complementary Items: *ComposiPro™ Brushes, VersaFlex™ Diamond-Coated Mylar Snap-On Discs and Strips, EP™ Flexible Pop-On Mylar Discs, VisionFlex® Perforated Diamond Strips, Gateway™ Solid Diamond Strips, ET® Composite Polishing Paste*

Handpiece Suggestions: *Forza™ ELXT Electric Motor, Forza™ M5 High Speed Electric Handpiece Attachment, Forza™ M1 Low Speed Electric Handpiece Attachment*

Recommended Speeds: *Polishers: 5,000-8,000 rpm; Final Polish: 2,000-4,000 rpm*

▶ Visit our website at **BrasselerUSA.com** or call **800.841.4522**.
In Canada call **800.363.3838**.

BRASSELER USA®

©2022 Brasseler USA. All rights reserved.