



BRASELER
USA®

TotalVac™

IRRIGATION SYSTEM



BY YOUR SIDE

THE CHALLENGE...


Multi-step, multi-product positive and negative pressure irrigation solutions, which are not able to be easily combined.

During a conventional root canal procedure, high-speed evacuation is used during the access preparation phase. The surgical suction is generally used during positive pressure irrigation. However, many clinicians prefer to use large volume negative pressure irrigation. Currently, the only available options are more expensive complete systems* with proprietary connections, needles and tips. All the current solutions require chair-side dental assistants, making it difficult for single-handed operation.

*Based on EndoVac™ Starter Kit retail pricing versus TotalVac Kit

THE SOLUTION...

TotalVac™ IRRIGATION SYSTEM

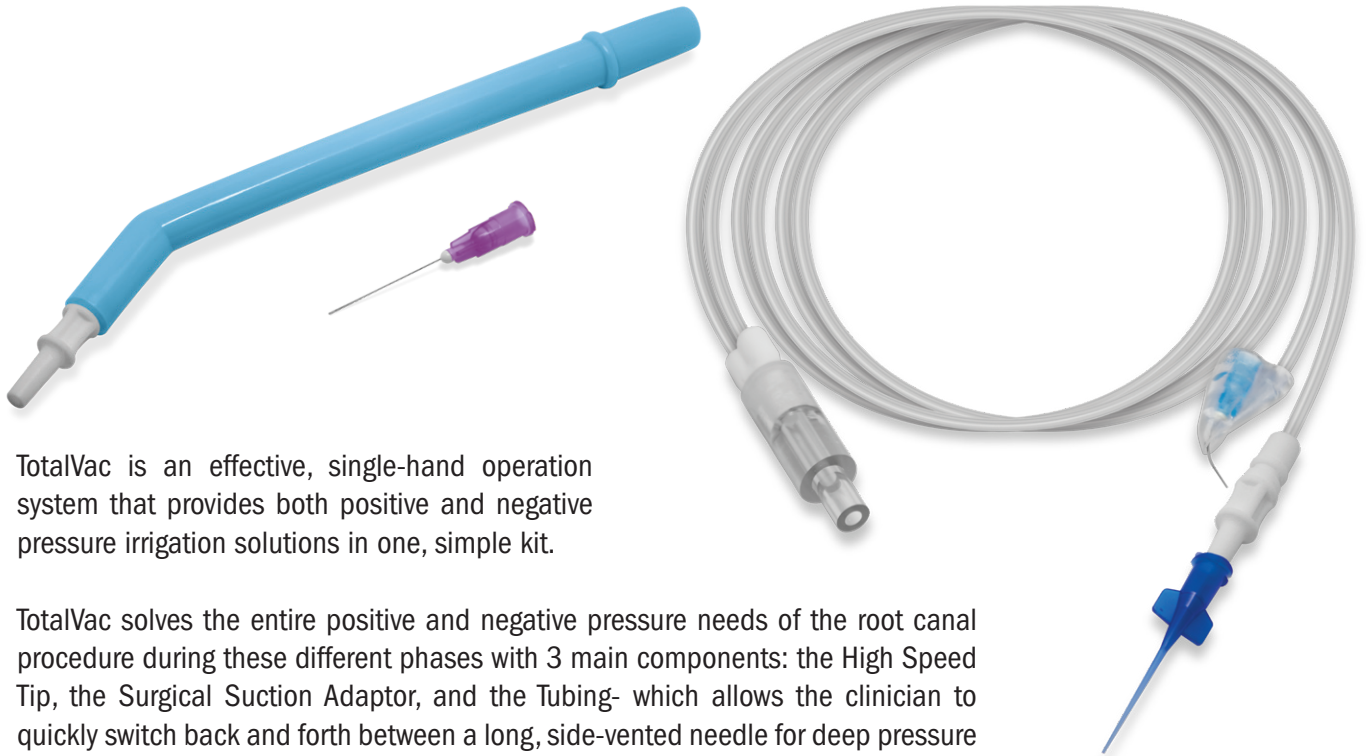

TotalVac is an effective, single-hand operation system that provides both positive and negative pressure irrigation solutions in one, simple kit.

TotalVac solves the entire positive and negative pressure needs of the root canal procedure during these different phases with 3 main components: the High Speed Tip, the Surgical Suction Adaptor, and the Tubing- which allows the clinician to quickly switch back and forth between a long, side-vented needle for deep pressure positive irrigation to a shorter needle on the device manifold, allowing for negative pressure irrigation while depositing irrigant inside the canal.

The TotalVac system allows the operator to address all operation suction needs, from access opening, to positive and negative pressure irrigation, all the way to drying the canal prior to obturation.

But more importantly, unlike other current systems on the market using their own proprietary connections and needle tips, TotalVac is an open-source system. It uses a smart, modular design allowing the operator to use their own preferred needles and designs.

FEATURES & BENEFITS


- ▶ Simple, convenient and easy to use
TotalVac is a completely modular, customizable system that allows you to easily perform both positive and negative pressure irrigation without the cost of purchasing expensive proprietary systems.
- ▶ The Complete Kit comes with all the items needed for clinical operation, but it's versatility allows it to be based on a clinician's own preferences and requirements
The standard luer lock attachment opens the kit up to possibilities for needle tip sizes, length and tip designs based on the operator preference. The included close-ended, side-vented needle provides safe and effective positive pressure irrigation as well as micro suction deep in the canal.
- ▶ TotalVac is the only system designed to allow for one-handed positive pressure irrigation
TotalVac is ideal for all operating environments but especially useful where access to a chair side dental assistant may be limited and one-handed positive pressure irrigation by the operator would be required.

PRODUCT OFFERING

TotalVac Kit


5027562U0

- 25 High Volume Evacuation and Surgical Suction Adaptor
- 25 TotalVac Tubing
- 25 Plastic Suction Tip
- 25 25 Gauge Short Delivery Tip
- 25 30 Gauge Side Vented Long Delivery Tip

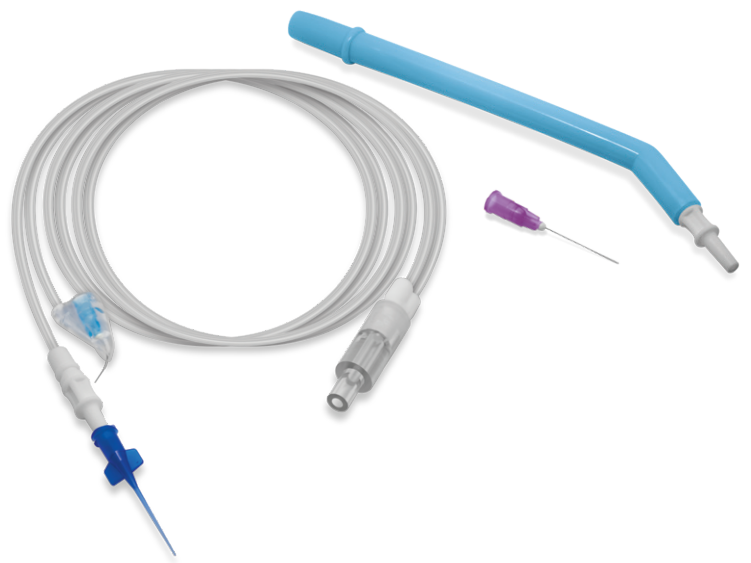


25GA Short Irrigation Tip (20/Pack)

5027569U0

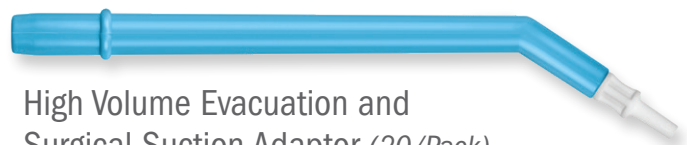

High Volume Evacuation and
Surgical Suction Adaptor (20/Pack)

5027566U0

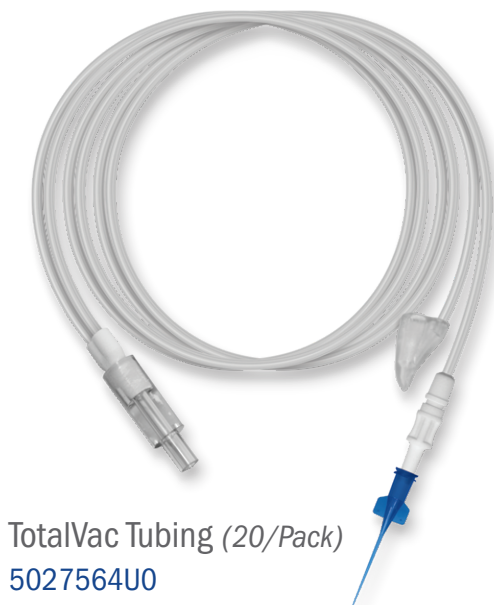


TotalVac Tubing (20/Pack)

5027564U0




Visit BrasselerUSADental.com/Discover-TotalVac
to read an article and watch an instructional video
presented by Dr. Allen Ali Nasseh, DDS, MMSc

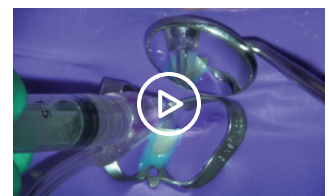



DIAMONDS

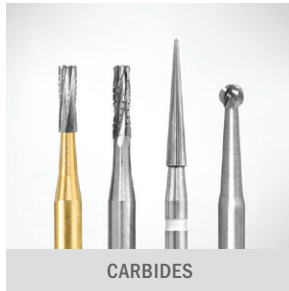

CARBIDES



POLISHERS



PROCEDURE SYSTEMS



BIOCERAMICS



Since 1976, our absolute focus has been to develop products and provide services that support the practice of core operatory procedures. Today, Brasseler USA® is the premier dental instrumentation company in North America offering the most extensive selection of dental and surgical instrumentation under one brand.



HANDPIECES



ENDODONTICS



SMALL EQUIPMENT



HYGIENE



HAND INSTRUMENTS



ORAL SURGERY



LABORATORY INSTRUMENTS



Visit our website at BrasselerUSA.com

To order call 800.841.4522 or fax 888.610.1937. In Canada call 800.363.3838

