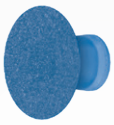
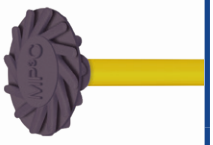
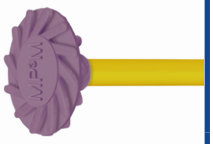
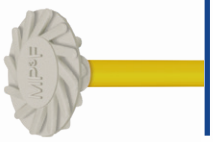
# Technique Guide: Removing Orthodontic Aligner Attachment & Bonding

## *Excess: Reduction & Removal*

| | Product Name | Order No. | Head Size ($^1/_{10}$ mm) | Cutting Length (mm) | Blade/Grit | Usage | Recommended Speed |
|---|---|---|---|---|---|---|---|
|  | Round-End Taper Trimming & Finishing Carbide | 7675 (H375R) | 016 | 8.0 | F 12-Blade Fine | Reduction and removal of excess adhesive with fine finishing carbide at a low speed with light hand pressure | 40,000 rpm |
|  | Long Flame Trimming & Finishing Carbide | H48L | 012 | 8.0 | F 12-Blade Fine | Reduction and removal of excess adhesive with fine finishing carbide at a low speed with light hand pressure | 40,000 rpm |
|  | Modified Round-End Taper Trimming & Finishing Carbide | H283K | 016 | 8.0 | F 12-Blade Fine | Reduction and removal of excess adhesive with fine finishing carbide at a low speed with light hand pressure | 40,000 rpm |
|  | VersaFlex™ Diamond-Coated Mylar Snap-On Disc | VFD1LG | 120 | – | Super-Coarse Green | Reduction and removal excess adhesive with a large mylar coarse disc at a low speed with light hand pressure | 7,000-10,000 rpm |
|  | EP™ Aluminum Oxide-Coated Mylar Pop-On Disc | EP1LR | 120 | – | Coarse Blue | Reduction and removal excess adhesive with a large mylar coarse disc at a low speed with light hand pressure | 10,000-12,000 rpm |

## *Residual: Removal & Polishing*

| | Product Name | Order No. | Head Size ($^1/_{10}$ mm) | Cutting Length (mm) | Grit | Usage | Recommended Speed |
|---|---|---|---|---|---|---|---|
|  | BrioShine™ Feather Lite™ Multi-Purpose Purple Polisher | MP3C | 100 | 1.6 | Coarse Dark Purple | Initial removal of residual adhesive on all surfaces, all restorative materials | 5,000-8,000 rpm |
|  | BrioShine™ Feather Lite™ Multi-Purpose Purple Polisher | MP3M | 100 | 1.6 | Medium Light Purple | Pre-polish and removal of final remnants on all surfaces, all restorative materials | 5,000-8,000 rpm |
|  | BrioShine™ Feather Lite™ Multi-Purpose Purple Polisher | MP3F | 100 | 1.6 | Fine Gray | Final polish, creates enamel-like finish on all surfaces, all restorative materials | 2,000-4,000 rpm |

*Complementary Items:* Forza V3 Piezoelectric Ultrasonic Scaler, EPM21 Pop-On Disc RA Mandrel, 310.21 Snap-On Disc RA Mandrel, K0315 VersaFlex Intro Kit (includes assortment of super coarse, coarse, medium, and fine discs in two sizes), EP200 EP Polishing Mylar System Kit (includes assortment of coarse, medium, fine and extra fine discs in two sizes)

*Handpiece Suggestions:* Forza Pro 1:1 Latch Electric Low Speed Handpiece Attachment
Forza Pro 1:5 FG Electric High Speed Handpiece Attachment

▶ Visit our website at **BrasselerUSA.com** or call **800.841.4522**.
In Canada call **800.363.3838**.

BRASSELER USA®

©2023 Brasseler USA. All rights reserved.