BRASSELER USA®

HIGH SPEED COUPLER

INSTRUCTIONS FOR USE

Read before use. Keep in a safe place for future reference.

CAUTION: Federal (US) law restricts this device to sale by or on the order of a licensed professional

Brasseler U.S.A. Dental, LLC hereinafter will be referred to as “Brasseler USA.”

Safety instructions are intended to avoid potential hazards that could result in personal injury or damage to the device. Safety instructions are classified as follows in accordance with the seriousness of the risk.

Class	Degree of Risk
⚠ **Warning**	Hazard that could result in serious injury or damage to the device if the safety instructions are not correctly followed
⚠ **Caution**	Hazard that could result in light or moderate injury or damage to the device if the safety instructions are not correctly followed.

Symbols	
SN	Serial Number
REF	Reference Number
QTY	Quantity
	Instructions for Use

⚠ Caution

- Do not modify or alter coupler.
- If handpiece leaks air or water and/or does not run up to speed, check o-rings or seals on coupler.
- Never attempt to remove the coupler while running handpiece.
- If any abnormalities are found, stop using immediately and contact Brasseler USA.

SPECIFICATIONS

ITEM CODE	DESCRIPTION	WEIGHT	LIGHT SOURCE
5027001U0	BR-4 Non-Optic 4-Hole Coupler	43 g	N/A
5027293U0	BR-5 5-Hole Coupler	43 g	Glass Rod Transfer
5027294U0	BR-6 LED 6-Pin Coupler	43 g	LED
5027002U0	BR-K5 5-Hole Coupler	26 g	Glass Rod Transfer
5027003U0	BR-K6 LED 6-Pin Coupler	25 g	LED

AIR PRESSURE

Set air pressure at handpiece connection to 38 - 42 psi.

⚠ Caution

- Exceeding recommended air pressure may result in premature cartridge failure.

OPERATION

Connecting Coupler (Fig. 1)

1. Align tubes and pins on coupler with holes in delivery system tubing.
2. Thread tubing nut to bottom of coupler.
3. Hold coupler and finger tighten

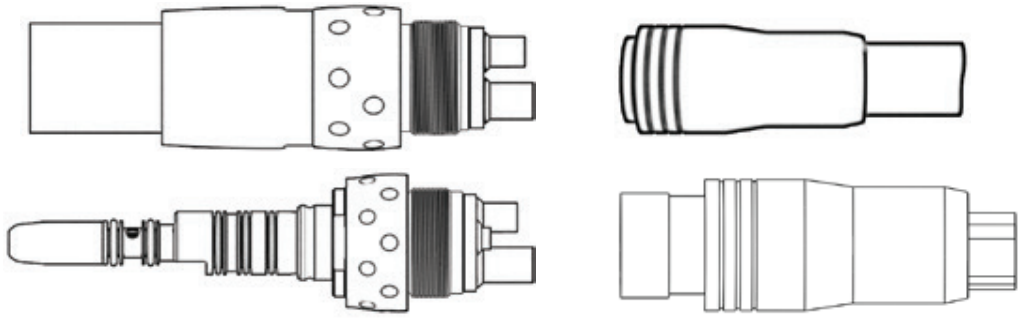
Fig. 1

Connecting Handpiece (Fig. 2)

1. Slide handpiece onto coupler.
2. Press firmly until it snaps into place.

Note: Handpiece should swivel freely with slight tension and sit flush on coupler.

Removing Handpiece (Fig. 2)

Brio 440 / Brio 440L: Pull back release ring on coupler and pull handpiece straight off.
Brio 440KL: Hold ring on coupler and pull handpiece straight off.

Note: Do not hold or pull on hose.

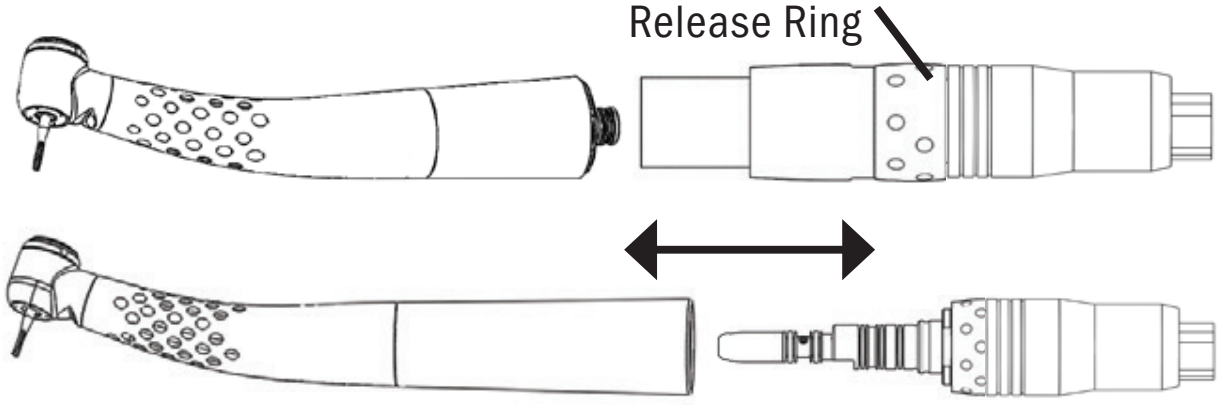
Fig. 2

MAINTENANCE

⚠ Caution

- Wear appropriate personal protective equipment during maintenance process. This should include eye protection and gloves.

Cleaning

1. Clear fiber-optic surfaces of particulates by spraying with air syringe.
2. Wipe fiber-optic surfaces with a soft cloth or cotton swab damp with alcohol.

Changing O-Rings
For BR-K5 and BR-K6 LED Couplers only.

⚠ Caution

- Improper care of o-rings may cause handpiece to malfunction.
- When changing/replacing coupler o-rings, always replace complete set of 5 pieces.

Note: Coupler o-rings may be lubricated with handpiece lubricant.

1. Use an explorer to gently remove old/worn o-rings.
2. Slide on new o-rings using only fingers to not cause damage.

Changing Coupler LED Bulb

Note: LED bulb is polarity sensitive. If bulb does not light after changing, repeat procedure and rotate bulb 180°.

For BR-K6 LED (Fig. 3)

1. Remove LED bulb cap by turning it counter-clockwise.
2. Remove bulb.
3. Install new bulb and replace cap, by finger tightening only.

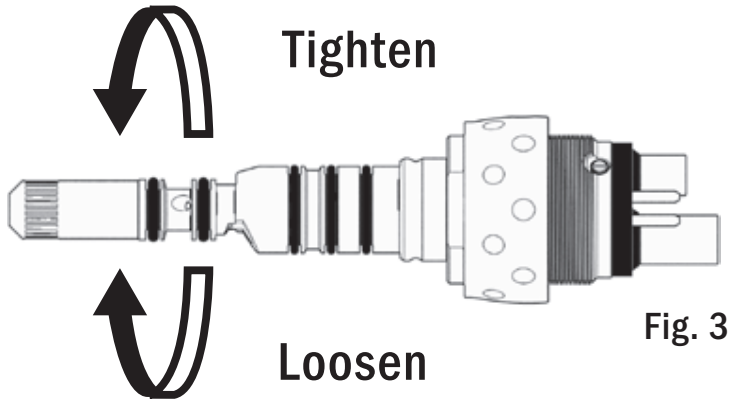
Fig. 3

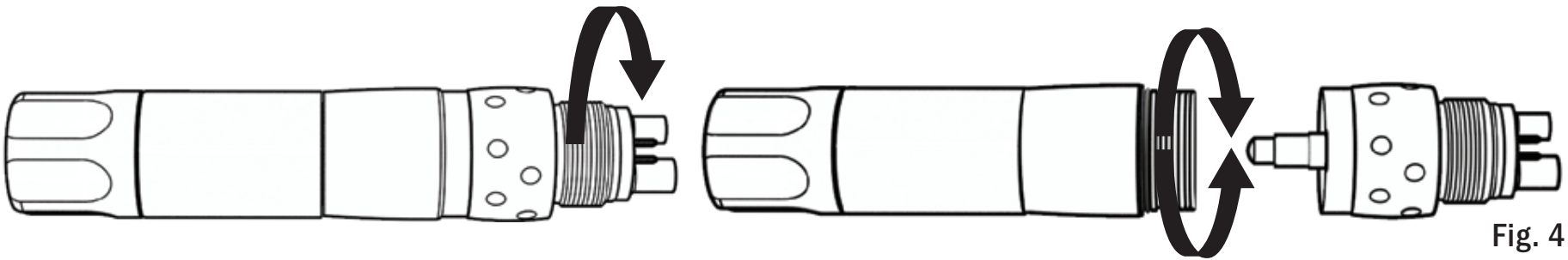
Fig. 4

For BR-6 LED (Fig. 4)

1. Use the coupling cap to remove coupling from the tubing.
2. While holding the coupling cap, rotate the coupling cap in the counter-clockwise of the arrow to loosen.
3. Remove the old light bulb from the hosing. Align and insert the connector of the new replacement bulb into the housing securely.
4. While holding the coupling cap, rotate the coupling cap in the opposite direction of the arrow to tighten.

REPLACEMENT PARTS & ACCESSORIES

5011053U0 BR-K4 Gasket

5012794U0 BR-6 LED / BR-K6 LED Gasket

5027125U0 BR-K6 LED Bulb

5011695U0 BR-K5 / BR-5 Gasket

5015061U0 BR-K5 / BR-K6 Coupler O-Ring Set

5027316U0 BR-6 LED Bulb

SERVICE

For items that require service, including warranty service, please send to Brasseler USA.

Brasseler USA
Attn: Dental Repairs
One Brasseler Blvd.
Savannah, GA 31419
800.841.4522 (in Canada, call 800.363.3838)

WARRANTY

These products are designed for use in a dental office and this warranty is not applicable to other uses. The products are warranted against manufacturing defects in materials and workmanship for 1 year from date of purchase.

Brasseler USA will repair or replace products at its sole discretion. This warranty shall be limited to replacement or repair of the products or their parts and shall not extend to any other claims, including but not limited to loss of profit, cost of removal or replacement, incidental or consequential damages or other similar claims arising from

use of this product. Damage to products resulting from acts of God, faulty installation, misuse, tampering, accident, abuse, negligence, alterations, unauthorized repairs or problems unrelated to materials and workmanship are not covered by this warranty.

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