BRASSELER USA®

FEATHER LITE™

KOIS UNIVERSAL POLISHERS

Developed with Dr. John Kois, a universal system designed for polishing all ceramic and composite restorations



K0430 Universal Polishing System
Simplifies integration and organization of all 6 polishers in one aluminum block

- Flexible unique shape; adapts to all tooth surfaces providing access to areas other shapes can't to create an enamel-like finish
- Fewer instruments and polishing steps; reduces chair time and inventory
- Reusable and autoclavable; cost-effective and sustainable
- Color-coded, diamond-impregnated polyurethane rubber, two-step polishing system for ease of identification and efficient polishing
- Preserves surface texture while polishing
- Metal latch mandrel for use with low speed handpieces*



Order No.	KUP1	KUP2	KUP3	KUP4	KUP5	KUP6
Head Size (1/10 mm)	140	140	140	55	55	55
Head Length (mm)	1.6	1.6	1.6	8.8	8.8	8.8
Grit	Medium	Medium	Fine	Medium	Medium	Fine
Color	Pink	Green	Gray	Pink	Green	Gray
Usage	Ceramic Pre-Polish	Composite Pre-Polish	Universal Final Polish	Ceramic Pre-Polish	Composite Pre-Polish	Universal Final Polish
21 - RA - Right-Angle Latch (22 mm)						
Qty/Pack	3	3	3	3	3	3

Indications:

- *Occlusal surfaces; pits and fissures*
- *Blends composite materials*
- *Removes lingual surface stains*
- *Lingual interproximal surfaces of anterior teeth and occlusal surfaces of posterior teeth*
- *Polishes ceramic after intra-oral occlusal adjustment*

Order online at **BrasselerUSA.com** or call **800.841.4522**.
In Canada call **800.363.3838**.

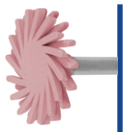
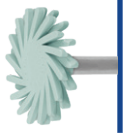
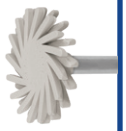
©2023 Brasseler USA. All rights reserved. Printed in the U.S.A.

B-5805-08.23

BRASSELER USA®

K0430 KOIS UNIVERSAL POLISHING SYSTEM

Developed with Dr. John Kois, Kois Center

	Instrument	Order No.	Head Size ($^{1}/_{10}$ mm)	Head Length (mm)	Grit	Usage
	Ceramic Feather Lite™ Polisher	KUP1	140	1.6	Medium Pink	Pre-polish all surfaces Occlusal surfaces; pits and fissures
	Composite Feather Lite™ Polisher	KUP2	140	1.6	Medium Green	Pre-polish all surfaces Occlusal surfaces; pits and fissures
	Universal Feather Lite™ Polisher	KUP3	140	1.6	Fine Gray	Final polish all surfaces Occlusal surfaces; pits and fissures Creates enamel-like finish to all restorations
	Ceramic Feather Lite™ Cup Polisher	KUP4	55	8.8	Medium Pink	Pre-polish cervical area of any teeth Lingual interproximal surfaces of anterior teeth
	Composite Feather Lite™ Cup Polisher	KUP5	55	8.8	Medium Green	Pre-polish cervical area of any teeth Lingual interproximal surfaces of anterior teeth
	Universal Feather Lite™ Cup Polisher	KUP6	55	8.8	Fine Gray	Final polish cervical area of any teeth Lingual interproximal surfaces of anterior teeth Creates enamel-like finish to all restorations

Complementary Items for Ceramic Adjustments:

- *8369DF and 8856DF Dialite® Diamonds*

Complementary Items for Composite Finishing:

- *8888.012 Fine Flame Diamond*
- *8392.016 Fine Mosquito Diamond*
- *ET6 or ET9 Needle T&F Carbides*
- *7408 (OS1) Football T&F Carbide*
- *8379.023 Fine Football Diamonds*
- *DET6F or DET9F Fine Needle Diamonds*
- *7901 Flame T&F Carbide*

****Brasseler Low Speed Handpiece Recommendations:***

- *Forza Pro 1:1 Latch Electric Attachment*

Brasseler High Speed Handpiece Recommendations:

- *Forza Pro 1:5 FG Electric Attachment*

Recommended speed: *4,000-6,000 rpm*
Final Polish: *2,000-3,000 rpm*

▶ Order online at **BrasselerUSA.com** or call **800.841.4522**.
In Canada call **800.363.3838**.

©2023 Brasseler USA. All rights reserved. Printed in the U.S.A.